

PRECISION  *MEDICAL*®

Oxygen Regulators



*The World's Most Compact and
Convenient Regulators!*

EASYREG OXYGEN REGULATOR

The World's Most Compact and Convenient Oxygen Regulators.



168708D



165415D



168715G

They're Compact and Lightweight...

The EasyDialReg is only 4.25 in. long and weighs a mere 7.8 oz. as a result of the offset gauge and multi-spring technology. The EasyGaugeReg is only 4.63 in. long and weighs 12 oz. with the contents gauge on the front and the flow gauge on the top.

They're User-Friendly...

Both types of regulators have easily accessible flow adjustment knobs making flow adjustment a cinch.

They're Safety is Built in...

The Easy Reg line of oxygen regulators feature a built-in pressure relief valve as well as a protected contents gauge. The low force Viton yoke seal provides a better tank-to-reg seal.

They're Safe...

The EasyDialReg has passed the ASTM G175 Ignition Sensitivity and Faulty Tolerance Test.

They're Accurate...

Unique design minimizes normal flow variation that occurs as tank pressure drops.


They're Flexible & Versatile...

The EasyDialReg comes in two (2) body designs yoke or nut style, and is available in five (5) flow ranges: 31-4000 cc/min, 1/32-4 lpm, .25-8 lpm, .50-15 lpm and .25-25 lpm. The EasyGaugeReg also comes in the two (2) body designs and two (2) flow ranges: 0-8 lpm and 0-15 lpm.

EasyDialReg Identification


DSDS	DSDR	DSR	NG	0706
168715DSDS shown Outlet DISS O2	168715DSDR shown Outlet DISS O2 90° pointing down	168715DSR shown Outlet Hose Barb 90° pointing down	168708DNG shown No Gauge	168715DSDS/0706 shown Extension used with DSDS for humidifier bottles

EASY DIAL REG
OXYGEN REGULATOR

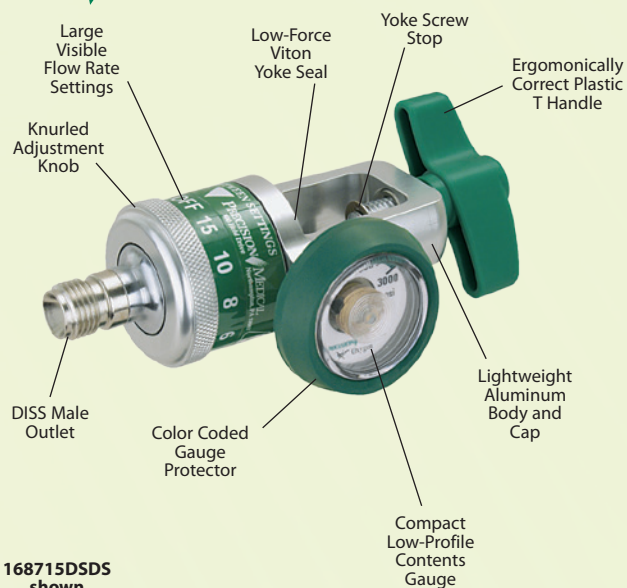

168715DSDS
shown

EASY GAUGE REG
OXYGEN REGULATOR

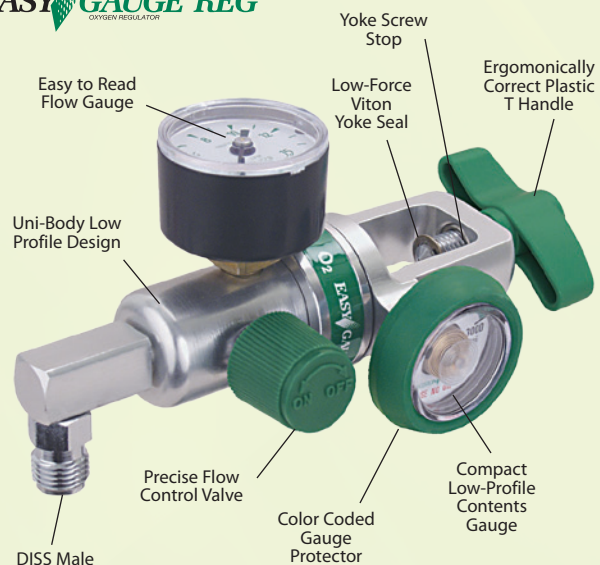

168715GE
shown

Ordering Information

EasyDialReg

Flow Range	Position Increments	CGA Connection	
		Yoke Style (870)	Nut Style (540)
0-4000 cc/min	OFF, 31, 62, 125, 250, 500, 750, 1000, 1500, 2000, 3000 and 4000	168704DCL	165404DCL
0-4 lpm	OFF, 1/32, 1/16, 1/8, 1/4, 1/2, 3/4, 1, 1-1/2, 2, 3 and 4	168704DL	165404DL
0-8 lpm	OFF, .25, .50, 1, 1.5, 2, 2.5, 3, 4, 5, 6 and 8	168708D	165408D
0-15 lpm	OFF, .50, 1, 1.5, 2, 3, 4, 5, 6, 8, 10 and 15	168715D	165415D
0-25 lpm	OFF, .25, .50, 1, 2, 3, 4, 6, 8, 10, 15 and 25	168725D	165425D

EasyGaugeReg

Flow Rate	Yoke Style (870)	Nut Style (540)
0-8 lpm	168708G	165408G
	168708GE	165408GE
0-15 lpm	168715G	165415G
	168715GE	165415GE

Specifications

Accuracy:	31 to 125 cc/min \pm 10cc .25 to 25 LPM \pm 10% of indicated flow
Transport/Storage Temperature Range:	-4°F (-20°C) to 140°F (60°C)
Inlet Pressure Range:	300 - 3000 psi (20.7 - 206.8 bar)
Ignition and Fault Tolerance:	meets ASTM G175-03
Warranty:	5 years

Specifications are subject to change without prior notice

Accessories & Replacement Parts

Yoke Seal

The low force Viton yoke seal provides a better tank-to-reg seal. It has a brass outer ring.



Ordering Information

1221-25 Yoke Seal (Pack of 25)

Disposable Tubing Nipple

Fits any Regulator with an Oxygen DISS fitting.



Ordering Information

0120 Yellow
0121 Black
0122 White
0124 Green

T-Handle

The T-Handle and Ring Retainer for the handle are parts for the Precision Medical Piston Regulators.



Ordering Information

502312 Ring, Retaining
504156 Handle Assembly, Yoke

Humidifier



PM500



PM510

Ordering Information

PM500 Disposable Humidifier
PM510 Autoclavable Humidifier
(50 per case)

Oxygen & Air Diaphragm Regulator 50 PSI Preset Regulators

These brass diaphragm regulators provide a minimum flow of 240 lpm at a constant pressure of 50 psi with DISS male outlet.



Ordering Information

Gas	CGA Connection	Style	Part Number
Oxygen	870	Yoke	870
	540	Nut	540
Medical Air	950	Yoke	504450
	346	Nut	504448

Color Gauge Guard Option

All models of Precision Medical's EasyReg can be ordered with any one of six different color gauge protectors, green being our standard. Different colors may be used to identify and segregate regulators by department. When ordering, please remember to add a dash (-) between the regulator part number and the desired color suffix. (Green gauge guard requires no suffix.) (Ex.: 168708D-Y)



PRECISION MEDICAL

Precision Medical, Inc.
300 Held Drive • Northampton, PA 18067
Phone: 610-262-6090 Toll-Free Phone: 800-272-7285
Fax: 610-262-6080 Toll-Free Fax: 800-352-1240
www.precisionmedical.com